

Plug in the USB mouse



Plug in power

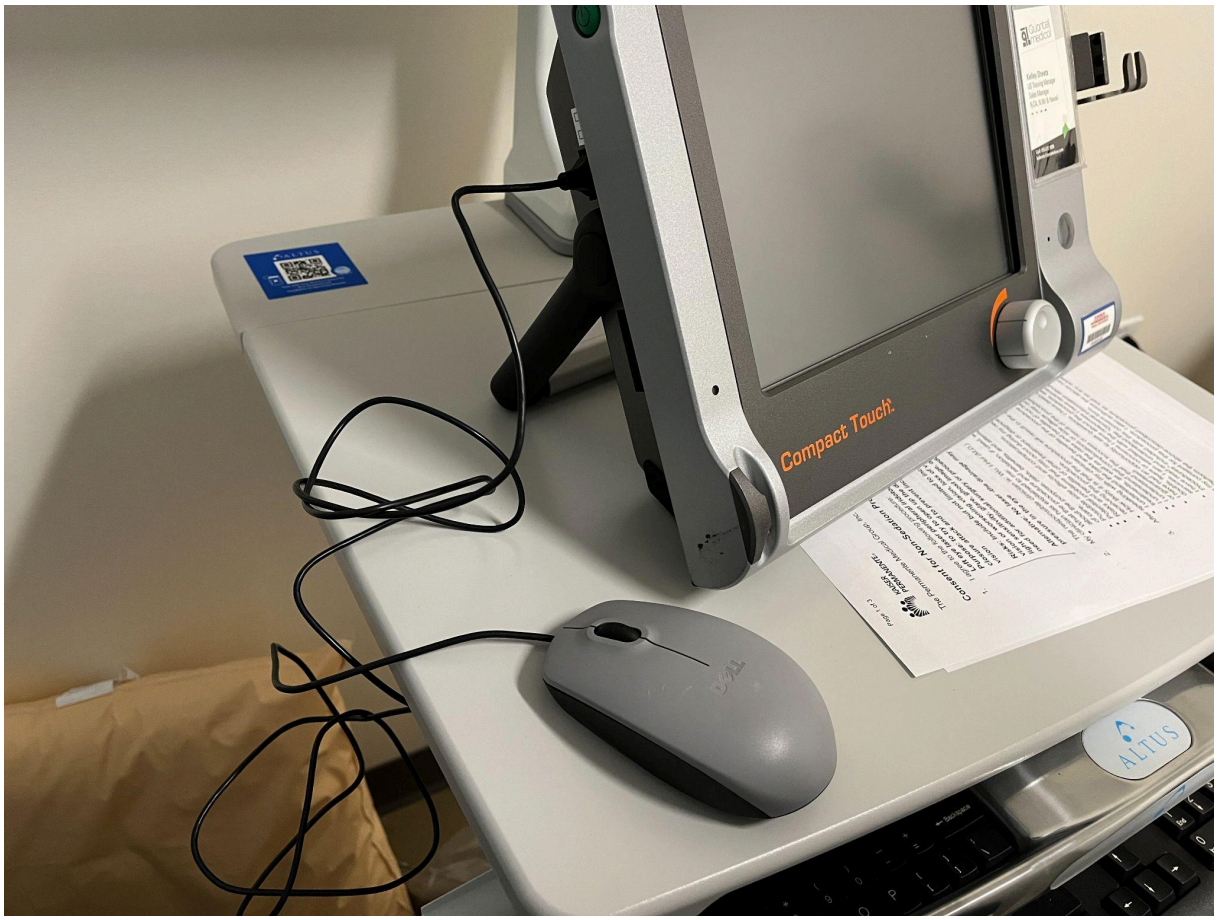

Connect the B-scan probe

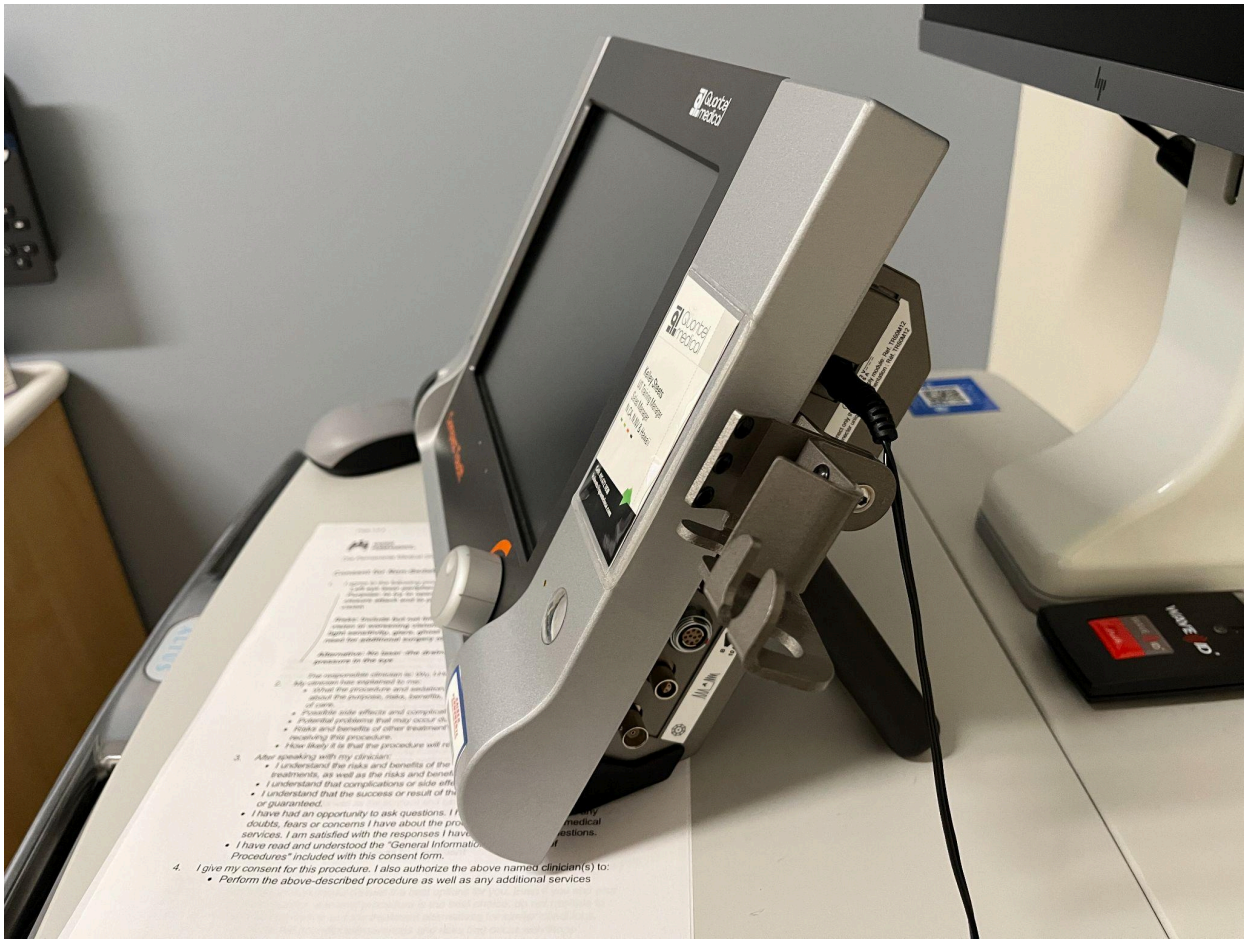

Turn on power with rocker switch under metal guard



Turn on device with green button




Press here to run B-scan



Click on footpedal icon to activate probe

Click on Gain if you want to adjust it with the dial.

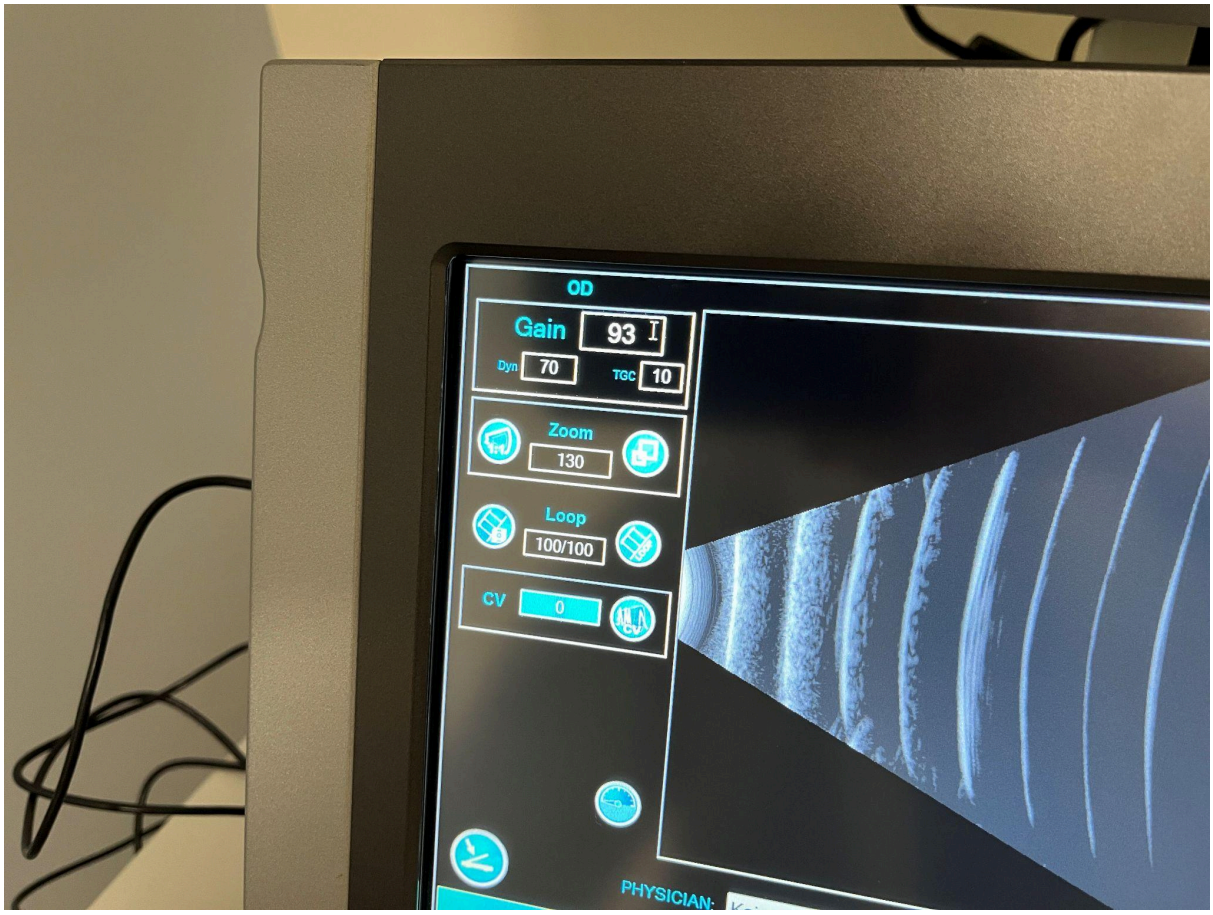
OD
Gain
93
Dyn
70
TGC
10
Zoom
130
Loop
100/100
CV
0
PHYSICIAN:


Kelley Sheets
US Training Manager
Sales Manager
N.CA, N.NV & Hawaii
Cell: 415.672.1838
ksheets@quantelusa.com
PHYSICIAN:
Kaiser Doc
Scan performed by:
<NONE>
Caliper
Marker
Area
Add to Report
pact Touch
KAISER PERMANENTE
Permanente Medical Group, Inc.