

Face W+WO IV Contrast Protocol

Revised: 11/11/21

Charge as: Face W+WO IV Contrast

Scanner: 1.5T or 3T


Coil: Head Coil/Flex Coil (flex coil may be needed if anterior head coil cannot clamp)

Plane	Weighting	Slice	Gap	FAT SAT	Scan Range
COR	T1	4mm	.4mm	None	Whole face/from tip of nose-4 th ventricle/ perpendicular to hard palate and perpendicular to nasal septum
AXIAL	T1	4mm	.4mm	None	Whole face/from glabella-larynx/parallel to hard palate and perpendicular to nasal septum
AXIAL	T2 FS	4mm	.4mm	Strong	Whole face/from glabella-larynx/parallel to hard palate and perpendicular to nasal septum
COR	STIR	4mm	.4mm	None	Whole face/from tip of nose-4th ventricle/ perpendicular to hard palate and perpendicular to nasal septum
SAG	T1	4mm	1mm	None	Whole face/from rt-lt pinna/parallel to nasal septum
AXIAL	Difussion	4mm	0mm	None	Whole face/from glabella-larynx/parallel to hard palate and perpendicular to nasal septum
Inject Contrast					
COR	T1 FS	4mm	.4mm	Strong	Whole face/from tip of nose-4th ventricle/ perpendicular to hard palate and perpendicular to nasal septum
AXIAL	T1 FS	4mm	.4mm	Strong	Whole face/from glabella-larynx/parallel to hard palate and perpendicular to nasal septum
SAG	T1 FS	4mm	1mm	Strong	Whole face/from rt-lt pinna/parallel to nasal septum

COR Face Image Planning:



AXIAL Face Image Planning:

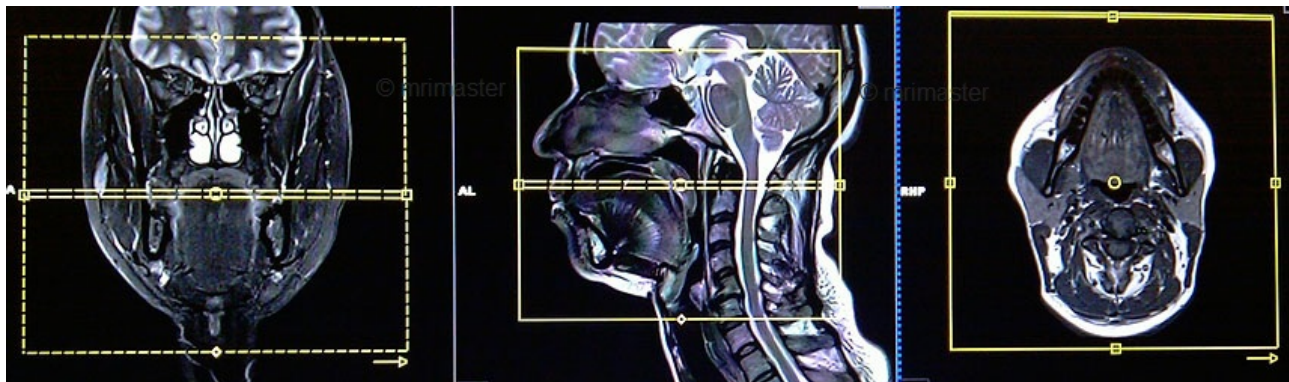

SAG Face Image Planning:

